


9/8/16

acl from previous day (containing unlabelled minipreps + ligations) was analyzed, found to have nothing of value.

• Parts contaminated?

Instead of ligating ~~ts~~ [pRlacI + ts purp (deg)] + [R011-clpXP CI]
we are now doing:

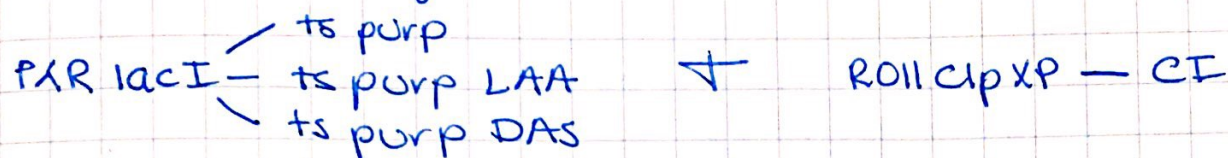

Digest (to make ends sticky)

1. ~~deg~~ digested pR lacI, ts purp (DAS/LAA), other groups digested backbones (1A3, 1C3, 1K3, 1T3) + R011-clpXP and CI

Calculations:

Parts were miniprepped at following concentration:

ts purp: ~~140.1~~ ^{140.1} ng/μL

ts purp DAS: 128.4 ng/μL

ts purp LAA: 130.0 ng/μL

pR lacI : 97.7 ng/μL

We need 1 ng (1000 ng) of DNA

ts purp:

$$\frac{0.1401 \mu\text{g}}{1 \mu\text{L}} = \frac{1 \text{ ng}}{x \mu\text{L}}$$

$$1/0.1401 = x$$

$$\boxed{7.138 \mu\text{L} = x}$$

ts purp DAS:

$$\frac{0.1284 \mu\text{g}}{1 \mu\text{L}} = \frac{1 \text{ ng}}{x \mu\text{L}}$$

$$1/0.1284 = x$$

$$\boxed{7.788 \mu\text{L} = x}$$

ts purp LAA:

$$\frac{0.1300 \mu\text{g}}{1 \mu\text{L}} = \frac{1 \text{ ng}}{x \mu\text{L}}$$

$$1/0.1300 = x$$

$$\boxed{7.69 \mu\text{L} = x}$$

pR lacI

$$\frac{0.0977 \mu\text{g}}{1 \mu\text{L}} = \frac{1 \text{ ng}}{x \mu\text{L}}$$

$$1/0.0977 = x$$

$$\boxed{10.235 \mu\text{L} = x}$$

continued to page

Signature

Date

Witnessed And Understood By

Date

PROPRIETARY
INFORMATION

8 Title

continued from page

amt of water needed -

$$\begin{aligned}
 \text{ts purp: } 20 - 7.138 (\text{DNA}) - 3 (\text{Enzymes}) &= 9.862 \text{ H}_2\text{O} \\
 \text{ts purp DAS: } 20 - 7.788 (\text{DNA}) - 3 (\text{Enzymes}) &= 9.212 \text{ } \mu\text{L H}_2\text{O} \\
 \text{ts purp LAA: } 20 - 7.69 (\text{DNA}) - 3 (\text{Enzymes}) &= 9.31 \text{ } \mu\text{L H}_2\text{O} \\
 \text{P} \lambda \text{R LacI: } 20 - 10.235 (\text{DNA}) - 3 (\text{Enzymes}) &= 6.765 \text{ } \mu\text{L H}_2\text{O}
 \end{aligned}$$

Protocol:

- 1) Do calculations to see how much μL of DNA you need to get ^{1 μg} μg
- 2) Determine which enzymes you need
 - Part one (P λ R LacI): E and S
 - Part two (ts purp (Q/PAS/LAA): X and P
- E and X are remixes, 2 μL of that are needed
- S and P are standard, 1 μL of it is needed
- 3) Put the ingredients into a microcentrifuge tube in this order: H₂O, DNA, enzymes. You should have 20 μL total.
- 4) Incubate at 37° for 1 hr, then at 80° for deactivation
 - ★ Use thermocycler + make diagram of tube placement because ink will wear off ★

8/9/16

- parts were ligated:

P λ R lacI ts purpP λ R lacI ts purp LAAP λ R lacI ts purp DAStransformations
for these were
started(10 μL of DNA was used)

R11 ClpXP C1 was miniprepped

• concentration: 53.7 $\mu\text{g}/\mu\text{L}$

• A260/A280: 1.86

9/12/16

transformation results:

- 3 phenotypes on each plate (P λ R lacI ts purp Q/LAA/DAS)
RFP, ts purp, normal

colony PCR on each phenotype

signature

Date

continued to page

continued from page

Purpose: to prepare biobricks for sequencingprotocol:

1) calculate amounts of DNA needed:

- There must be 12 μL total - 0.8 μL Primer, X DNA + rest
- 1 μg (1000 ng) DNA, so calculate w/ miniprep conc.

PAR lacI : 90.9 ng/ μL

$$\frac{90.9 \text{ ng}}{1 \mu\text{L}} = \frac{1000 \text{ ng}}{X \mu\text{L}}$$

$$11.0 \mu\text{L} = X$$

ts purple ①

$$\frac{77.4 \text{ ng}}{1 \mu\text{L}} = \frac{1000 \text{ ng}}{1 \mu\text{L}}$$

$$12.9 \mu\text{L}$$

ts purple ②

$$\frac{364.3 \text{ ng}}{1 \mu\text{L}} = \frac{1000 \text{ ng}}{X \mu\text{L}}$$

$$2.74 \mu\text{L} = X$$

ts purp ③

$$\frac{120.0 \text{ ng}}{1 \mu\text{L}} = \frac{1000 \text{ ng}}{1 \mu\text{L}}$$

$$8.33 = X$$

* total rxns:PAR lacI : 11.0 μL DNA + 0.8 μL primer + 0.2 μL H₂O

ts purple ①:

ts purple ②: 2.74 μL DNA + 0.8 μL primer + ~~12.9~~ 8.5 μL H₂Ots purple ③: 8.33 μL DNA + 0.8 μL primer + 2.9 μL H₂O

- 2) puncture rubber top with pipet tip. Put ingredients inside
- 3) Mail using prepaid envelope.

continued to page

continued from page

9/22

transformation of

- DAS ligation w/ 1c3
 - LAA ligation w/ 1c3
- } onto cam plates

- 1) Thaw materials on ice for 5 minutes
- 2) Put 10 μ L of ligation onto 100 μ L (1 tube) of cells
- 3) Flick the tube to mix
- 4) Put on ice for 30 minutes
- 5) Add 200 μ L of LB media / SOC
- 6) Incubate at 1 hr at 37°C
- 7) Plate 150 μ L of cells onto a plate, to a lawn
★ MAKE SURE YOU USE PROPER ANTIBIOTIC

- Used C2987, NOT NEB Beta cells

9/23

★ successful 😊

continued from page

10/4 : sequencing
EGFP-3p.F primer (10 μ M)

~~tube 1 tube 2~~
~~VF2-VF: RII ClpXP CI (nd first)~~
#3 EGFP : PAR GFP LAA
#4 EGFP PAR GFP DAS
#5 EGFP PAR GFP Φ

Final Plans:

- 1) PAR GFP Φ ClpXP CI (48.7 ng/ μ L)
- 2) PAR GFP DAS ClpXP CI (297.3 ng/ μ L)
- 3) PAR GFP LAA ClpXP CI (285.6 ng/ μ L)

calculations:

$$1) \frac{48.7 \text{ ng}}{1 \mu\text{L}} = \frac{1000 \text{ ng}}{x \mu\text{L}}$$

$$20.5 \mu\text{L} = x$$

no tag

$$2) \frac{297.3 \text{ ng}}{1 \mu\text{L}} = \frac{1000 \text{ ng}}{x \mu\text{L}}$$

$$3.363 = x \mu\text{L}$$

DAS

$$1.837 \text{ H}_2\text{O}$$

$$3) \frac{285.6 \text{ ng}}{1 \mu\text{L}} = \frac{1000 \text{ ng}}{x \mu\text{L}}$$

$$3.5 \mu\text{L} = x$$

LAA

$$1.7 \text{ H}_2\text{O}$$

tube AHA756 : PAR GFP LAA ClpXP CI
tube AHA757 : PAR GFP DAS ClpXP CI
tube AHA753 : PAR GFP Φ ClpXP CI

Backbone Switch for:

- 1) PAR LacI (PAR L)
- 2) RODII ClpXP CI (RCpCI)
- 3) PAR LacI GFP LAA RODII ClpXP CI (pLAA CpCI)
- 4) PAR LacI GFP DAS RODII ClpXP CI (pDAS CpCI)

1) Digest

(not digesting RCpCI, not enough)

Amounts:

PAR LacI: 2 μ L E + 1 μ L P + 11.0 μ L DNA + 6 μ L H₂O = 20 μ L

PLAACpCI: 2 μ L E + 1 μ L P + 17.0 μ L DNA = 20 μ L

PDASCpCI: 2 μ L E + 1 μ L P + 17.0 μ L DNA = 20 μ L

conc. 90.9



Order: Water:

DNA

Enzymes

10/14:

redigested the above mixtures, due to inconsistent levels of 10/13 miniprep

Ligation calculation

Parts:

PAR Lin 3TS) 4a+10GS (3TS: 3200)

PLAACpCI (in AT3) above + ~~753~~ 753 + 2974

pDASCpCI (in AT3) above + 153 + 2974

(AT3: 3450)

Part	length of F	length of V	total		μ L
Vector 13	20	2100	2100	4.7	1
PAR LacI	1114	3200	4314	0.23	(2100) $\times 3 = 6.12$
PLAA	829 484	3450	829	0.12	(3.91) $\times 3 = 11.75$
pDAS	829 484	3450	<u>829</u>	0.12	" " = 11.75

10/19

transformation of ligation (see above) onto CAM plates
sequencing

ROBII CIPXP CI in 1AK3, ~~10 beta cells~~ : 175.8 $\frac{\text{ng}}{\mu\text{L}}$ (1.85)
PAR LacI cFP & CIPXP CI in 1C3 : 79.1 $\frac{\text{ng}}{\mu\text{L}}$ (1.69)

Calculations:

ROBII CIPXP CI:

$$\frac{175.8 \text{ ng}}{1 \mu\text{L}} = \frac{1000 \text{ ng}}{x \mu\text{L}}$$

5.68 μL

PAR...XP CI:

$$\frac{79.1 \text{ ng}}{1 \mu\text{L}} = \frac{1000 \text{ ng}}{x \mu\text{L}}$$

12.64 μL

ROBII: 5.68 μL DNA + .8 prim. + 5.52 H₂O
PAR: 12.64 μL DNA + .8 prim

ROBII CIPXP CI w/ VF2 : ADF954
" w/ VR : ADF956

PAR... CI w/ VF2 : ADF953
" w/ cGFP : ADF955
w/ VR : AHA758

10/16

PAR LacI grew, pLADAS CIPXP CI + pLGLAACIPXP CI did not.
We are re-digesting + ligating them.

continued to

Signature	Date
Witnessed And Understood By	Date